

MULTIPOL ADVANCED PLUS

Automatic Grinder - Polisher

Multipol Advanced Plus is an automatic grinding and polishing device used in metallographic sample preparation. Thanks to its servo technology, it allows you to achieve excellent results even in the most difficult samples.

- | Single or dual disc design
- | Easy cleaning thanks to extendable and flow-adjustable taps
- | Optional dosing system
- | The sample holder unit can be used for sensitive coatings that need to be prepared on one side while remaining stable during operation.
- | Self cleaning basin system



/ SPECIFICATIONS

50-650 RPM

Rotation Speed


0-200 RPM

Head Speed

10-60 N

Pressure Force

/ PRODUCT DETAILS

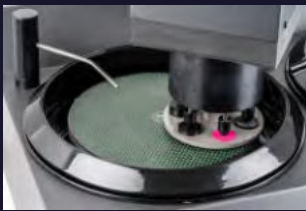

/ TECHNICAL SPECIFICATIONS

Disc Options

Ø200/250/300mm

Sample Capacity

6xØ40-25mm
3xØ50mm

Control

Color touch panel

Rotation direction

Both direction

Motor (disc)

AC/0.75kW

Motor (head)

Servo/0.75kW

Dimensions (w x d x h)

585x695x700mm *
760x695x700mm **

Weight

95kg *
120kg **

Power Supply

200-240V
50/60Hz

*Single Disc ————— **Double Disc